

The Procedure of Installation :

The first step: overview the product (bluetooth device) + (conductive mask sheet).

Step 2: Using USB cable to charge the smart device for the first use, red indicator light when charging, blue indicator light when fully charged (Note that device cannot be used during charging).

Step 3: When the battery is fully charged and the blue indicator light is on, please using mobile phone to download app according to the below of procedure;

iPhone : Please visit APPLE STORE downloads the “**smart beauty**” APP and install it.

Android Phone : Please scan QR code as below and install it.

Turn on the "Bluetooth" and "Location" of the mobile phone before the APP download, the place to be ticked when downloading; please tick to” allow “.



The procedure of APP operation as follows.

(1). After the APP is installed, the phone screen displays "Smart Beauty".

(2). The bluetooth device is powered on and turned on (the blue light flashes once and again), and it is recommended to fully charge it for the first time use.

(2). The Bluetooth of the mobile phone should be turned on.

(3). Open the smart beauty APP,

A. Click "Add Device" as shown in Figure 1;



Figure 1

B. Choose to add “smart mask, eye mask, shoulder and neck membrane, pleura, peritoneum, waist membrane, private membrane...etc” , depend on what product is added and click on the corresponding nam as shown in Figure 2;

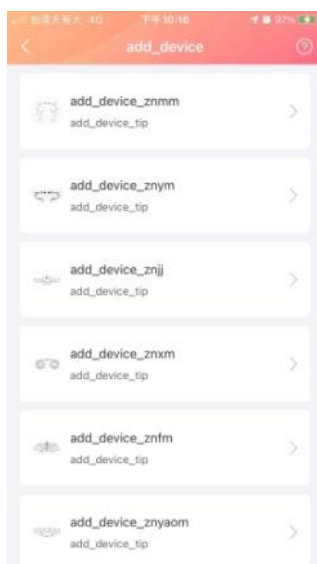


Figure 2

C. Click to search as shown in Figure 3

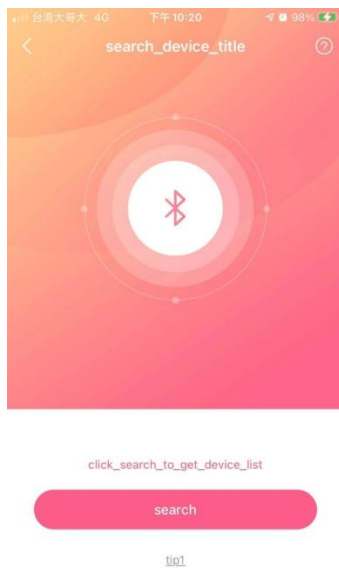


Figure 3

D. Click on the connection as shown in Figure 4.

Enter Figure 5 and enjoy the 12-minute experience.



Figure 4

E. Mode Selection :

According to the user's experience to choose the right mode for therapy, each time use

for 12 minutes , please refer to Figur 5.

If change mode , the current intensity will be returned to “1” level , user need to adjust the intensity properly again.

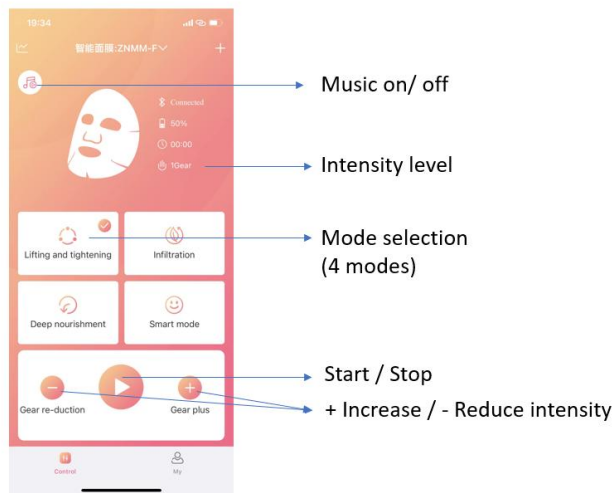


Figure 5

F. Figure 6 and Figure 7 are the trend graph indicators after the 12-minutes of experience (for reference only);

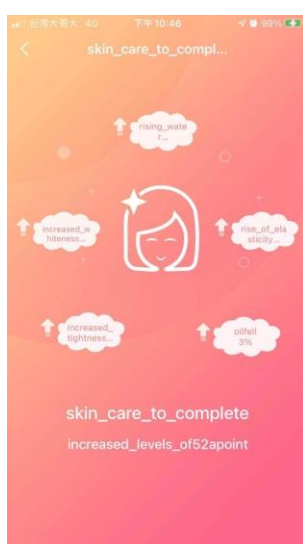


Figure 6

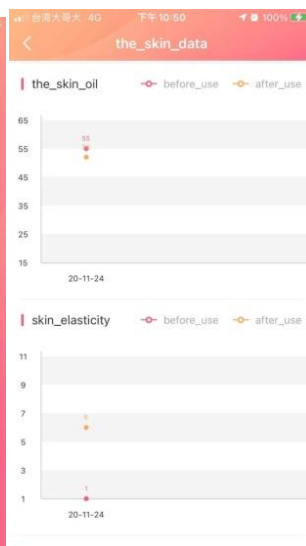


Figure 7

Notices:

1. Bluetooth connection is not available

(1) After the mobile phone downloads the APP, open the APP, a dialog box will pop up, select "Always allow" to enable bluetooth, if you don't click Allow, please choose to turn on Bluetooth through the phone settings.

(2) For the first time of use, please read the manual to charge the device. The red light will be on when charging, and the blue light will be on when fully charged.

It is expected to charge for 2 hours; note that the charger voltage is 5V/2A or less.

(3) Turn on the device, the blue light will flash, and the opened APP will automatically connect and click to connect. If the connection is successful, the APP will be displayed (there is a MAC number on the mobile phone APP), if the connection is not successful, please check as follows

A. First at all ,check whether the bluetooth device is turned on (whether there is a blue light flashing) and confirm whether the device has power and the phone's Bluetooth is turned on.

B. Whether the app is open or not.

C. Switch phone to try again.

D. If the device stop suddenly during the experience , it is might be caused by the following reasons, please handle it as follows:

a. Please turn off the host and close the APP.

b.To confirm whether two mobile phones have opened the APP at the same time. It is recommended to connect one by one. Start the APP and turn on the device.

Start using after connection.

2. Pulse Wave

(1) After the bluetooth device is connected to the APP, please attach device to the two contacts of the conductive mask sheet , adjust the level of massage intensity and enjoy a 12-minute experience . Note that the intensity is adjusted from low to high (1-16 levels).


(2) The intensity of the pulse wave can be adjusted according to your own feelings.

everyone's tolerance of bioelectric microcurrent stimulation is different, so there may be very few people who have a bioelectric stimulation reaction, which is normal; if you feel obvious discomfort during use, you should immediately lower the intensity or stop using it.

3. Data display

(1) The smart device combines with the conductive mask sheet to complete a perfect experience.

(2) At the end of 12 minutes, there will POP UP a picture of skin changes immediately and record a trend chart report, you can refer to the trend report, After you use it for the second time, the graph of curve will show data in skin care report.

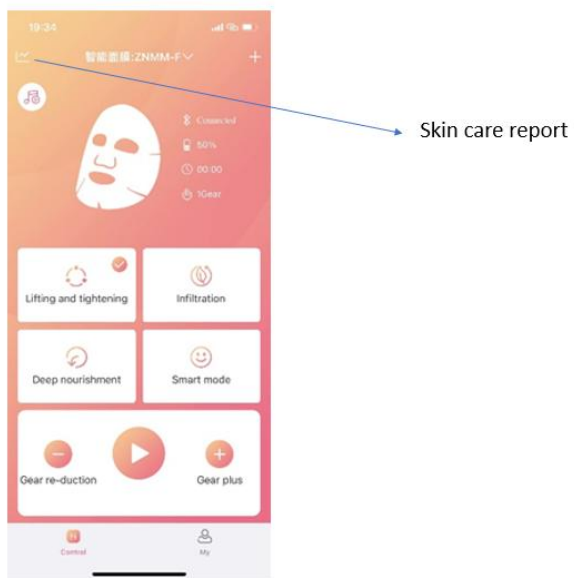

4. Maintenance

- (1) Make sure to clean the magnetic contacts of bluetooth device after use.
- (2) When the bluetooth device is in use, it must not fall on the ground (especially on concrete floors and tiles).

5. How to use conductive mask sheet?

Figure A. Put the side with two **white insulating points** toward the skin.



Figure A

Figure B

Figure B.C.D. The bluetooth device is attached on **grey color of two magnetic**

contacts on the silver paste surface.



Figure



Figure D

Precautions:

1. In the process of using the smart facial mask, you can feel the slight "electric light" on the two outer sides of the face and eyes with a slight dazzling sensation when you "close your eyes", which are all normal reactions.
2. Everyone has different tolerance to bioelectrical microcurrent stimulation, so there may be very few people who have a bioelectrical stimulation reaction. This is normal. If you feel uncomfortable during use, you should immediately lower the intensity of micro current or stop using it. ;
3. People with light-sensitive or allergic skin should use it with caution (preliminary judgment of light-sensitive: in summer, the skin will be red, the skin stratum corneum is particularly thin, there are red blood streaks, there will be a burning sensation after a few minutes in the sun, and the skin will change to very dark, basically light sensitive);
4. For those with sensitive skin, please use the essence to do a local skin sensitivity test before use. Use it without irritation. If you have any discomfort or persistent redness, you should stop using it immediately;

5. Pregnant women, patients with high blood pressure, heart disease, face metal implants, heart stents, organic diseases and skin wounds are prohibited;

6. Please put it out of the reach of infants, do not put it in a place where it is exposed to sunlight or high temperature;

7. After each use, wipe the remaining essence on the smart device to keep the bluetooth device clean;

A. When the bluetooth device is in use, it must not fall on the ground (especially on concrete floors and tiles);

B. Ensure that the contacts are wiped clean after use, and placed in a dry place; to avoid rusting of the magnet.